

# BIOSENSOR ASSAY ENHANCEMENT THROUGH AC ELECTROKINETIC MICROSTIRRING

M. Sigurdson, H. Feldman, and C. D. Meinhardt\*  
University of California, Santa Barbara  
Santa Barbara, CA 93106

## ABSTRACT

A microstirring tool is developed to improve binding rates in biosensor devices by an order of magnitude. AC electric fields generate fluid motion through the well documented but unexploited phenomenon, Electrothermal Flow, where the circulating flow redirects or stirs the fluid, providing more binding opportunities between suspended and wall-immobilized molecules. Simulations and experiments show improvement by a factor of 6-13. This is applicable to diffusion limited assays of different formats, from handheld microfluidic to benchtop microarray. Electrothermal microstirring is a tool that can be incorporated in existing devices, to speed up, for example, food-borne pathogen detection, or portable medical diagnostics.

## 1. INTRODUCTION

### 1.1 Immunoassays

Immunoassays, which rely on specific antigen-antibody binding for identification of proteins in a sample, have applications in both clinical laboratories for medical diagnostics and treatment monitoring, and in research laboratories for highly multiplexed testing, such as for biomarker identification. Additionally, portable immunoassays can be used for detection of environmental or food-borne pathogens, or soldier medical diagnostics. In these cases, throughput is a key consideration. One factor that can limit test duration is diffusion rate of analyte to the reporter. This is particularly true for high sensitivity ELISA tests (Stenberg; Nygren 1988). An incubation step of minutes to hours allows diffusion-limited binding to reach detectable levels.

In response to the needs for increasing throughput, portability, and sensitivity, new formats for miniaturized immunoassays have developed dramatically in recent years. These include spotted microarrays, common for “gene chips”, and now used for “protein chips” (e.g. the ProtoArray from Invitrogen, Carlsbad, CA), and various forms of Lab-on-a-chip devices which can perform fluid processing and detection steps on a single chip (Northrup; Jensen et al. 2003). Small length scales permit small sample sizes and shorter test times; on-chip sample preparation reduces fluid handling steps. Though greatly aided by their small length scales, these assays

can still be limited in response by diffusion of the analyte to an immobilized ligand. In this paper, we present a method where AC electrokinetics is used to enhance the performance of heterogeneous assays. Here electrothermal forces are used to micro-stir the analyte near a functionalized surface, increasing the rate of transport to the surface. This addresses the need for faster assays by offering a tool that is adaptable to a wide variety of assay configurations and reduces the incubation time required to perform a specific test while maintaining its sensitivity.

### 1.2 AC Electrokinetics

AC electric fields in a microchannel or microcavity can generate forces on both the fluid itself, through AC electroosmosis and electrothermal forces, and on suspended particles through dielectrophoresis. These forces are particularly important at the micron scale where buoyancy loses its influence and where high field gradients are possible with relatively low voltages. Electrothermal flow arises from Joule heating, which produces temperature inhomogeneities and therefore temperature-dependent conductivity and permittivity inhomogeneities, which interact with the applied electric field to produce fluid motion. The electrothermal body force is given by

$$\vec{F}_{ET} = -0.5 \left[ \left( \frac{\nabla \mathbf{s}}{\mathbf{s}} - \frac{\nabla \mathbf{e}}{\mathbf{e}} \right) \cdot \vec{E} \frac{\mathbf{e} \cdot \vec{E}}{1 + (\mathbf{w} \mathbf{t})^2} + 0.5 |\vec{E}|^2 \nabla \mathbf{e} \right] \quad (1)$$

(Ramos; Morgan et al. 1998), where  $\mathbf{t} = \mathbf{e}/\mathbf{s}$  is the charge relaxation time of the fluid medium and the incremental temperature-dependent changes are


$$\nabla \mathbf{e} = \left( \frac{\partial \mathbf{e}}{\partial T} \right) \nabla T, \quad \nabla \mathbf{s} = \left( \frac{\partial \mathbf{s}}{\partial T} \right) \nabla T. \quad (2)$$

Electrothermal forces are particularly useful in highly conductive (e.g. saline buffer) solutions, where increased fluid conductivity yields more heating of the fluid at lower applied voltage.


Report Documentation Page				Form Approved OMB No. 0704-0188	
Public reporting burden for the collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington VA 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to a penalty for failing to comply with a collection of information if it does not display a currently valid OMB control number.					
1. REPORT DATE <b>01 NOV 2006</b>		2. REPORT TYPE <b>N/A</b>		3. DATES COVERED <b>-</b>	
4. TITLE AND SUBTITLE <b>Biosensor Assay Enhancement Through Ac Electrokinetic Microstirring</b>				5a. CONTRACT NUMBER	
				5b. GRANT NUMBER	
				5c. PROGRAM ELEMENT NUMBER	
6. AUTHOR(S)				5d. PROJECT NUMBER	
				5e. TASK NUMBER	
				5f. WORK UNIT NUMBER	
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) <b>University of California, Santa Barbara Santa Barbara, CA 93106</b>				8. PERFORMING ORGANIZATION REPORT NUMBER	
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES)				10. SPONSOR/MONITOR'S ACRONYM(S)	
				11. SPONSOR/MONITOR'S REPORT NUMBER(S)	
12. DISTRIBUTION/AVAILABILITY STATEMENT <b>Approved for public release, distribution unlimited</b>					
13. SUPPLEMENTARY NOTES <b>See also ADM002075., The original document contains color images.</b>					
14. ABSTRACT					
15. SUBJECT TERMS					
16. SECURITY CLASSIFICATION OF:			17. LIMITATION OF ABSTRACT <b>UU</b>	18. NUMBER OF PAGES <b>4</b>	19a. NAME OF RESPONSIBLE PERSON
a. REPORT <b>unclassified</b>	b. ABSTRACT <b>unclassified</b>	c. THIS PAGE <b>unclassified</b>			

## 2. ELECTROTHERMAL VELOCITY MEASUREMENTS

In order to measure electrothermal velocity, a planar flow device was fabricated (Fig 1). A large flow cavity was through-etched in a thin (200  $\mu\text{m}$ ) silicon wafer; sidewall electrodes were deposited with a gap between them (around 90  $\mu\text{m}$ .) The silicon device was sandwiched between glass and filled with fluorescent tracer particles in 0.05M KCl. Micro-PIV was used to measure particle velocity and then estimate fluid velocity (Wang; Sigurdson et al. 2005). The measured velocity (Fig. 2) is on the order of 100  $\mu\text{m/s}$  for 6  $V_{\text{rms}}$  applied to the electrodes.

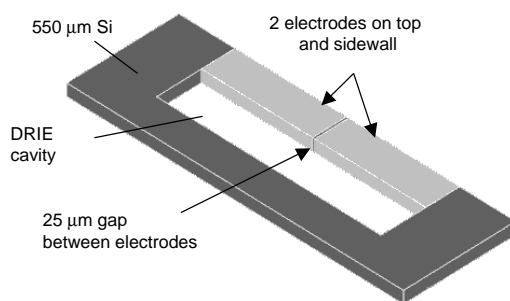


**Fig. 1** Test device used to measure 2D flow created by an inhomogeneous AC electric field. The device is sandwiched between glass, and the cavity is filled with 0.05 M KCl solution seeded with fluorescent polystyrene particles. An AC driving signal produces an electric field in the fluid which generates fluid motion.

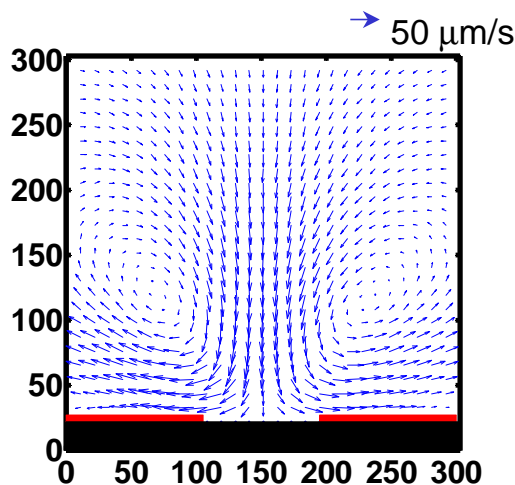

**Fig. 2** Fluid velocity field measured in the device shown in Fig. 1. Images of the fluorescent particles moving in the fluid are captured at discrete times. Micro-Particle Image Velocimetry is used to determine particle velocity field. Data from two different size particle velocity fields are combined to estimate the fluid velocity. The electrodes are driven at 7  $V_{\text{rms}}$  and 200 kHz. The AC field produces two counter-rotating vortices of characteristic velocity 100  $\mu\text{m/s}$ .

## 3. FLOW & BINDING SIMULATIONS


The finite element package *Femlab* (Comsol; Stockholm, Sweden) was used to solve the electric and temperature fields, and the electrothermal force (Eq. 1). The result is used as a forcing term in the Navier Stokes equations, and the velocity field is solved (Sigurdson; Wang et al. 2005). The simulation results are shown in Fig. 3, and compare well with experimentally measured velocity (Fig. 2). In order to describe how this circulating fluid motion can enhance binding rate of an immunoassay, we solved for the convection and diffusion of a scalar in the microchannel, subject to a binding boundary condition. In the passive case (Fig. 4a), a depleted cloud surrounds the binding region, and the binding rate is reduced. When electrothermal microstirring is turned on, (Fig. 4b) the depleted concentration is circulated away, and fresh analyte is circulated past the sensor. This increases binding rates.

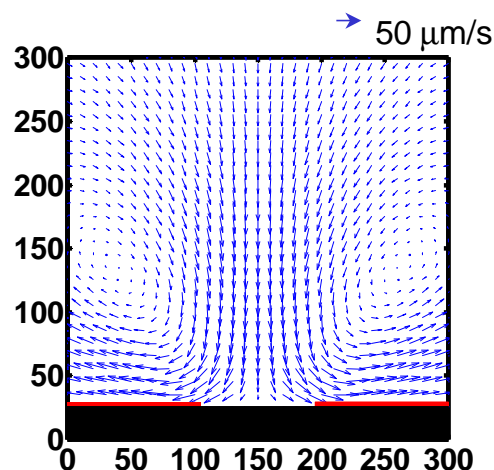


**Fig. 3** Simulation of electrothermally-driven fluid motion with an applied voltage of 7  $V_{\text{rms}}$ . The velocity of the electrothermally-driven flow is of order 100  $\mu\text{m/s}$  and is characterized by a pair of counter rotating vortices. This flow closely matches the velocity field that was experimentally measured using micro-PIV (Fig. 2).

## 4. BINDING EXPERIMENTS


To test electrothermal microstirring with a heterogeneous streptavidin-biotin reaction, a different geometry was used. Two electrodes separated by 55  $\mu\text{m}$  were deposited on silicon with an  $\text{SiO}_2$  isolation layer (Fig 5). Because the electrothermal flow is driven by a temperature *gradient* (Eq. 2), and the Joule heating provides a warm spot within the fluid, it is important to ensure sufficient heat removal at one or more walls in order to maintain a temperature gradient. In this device

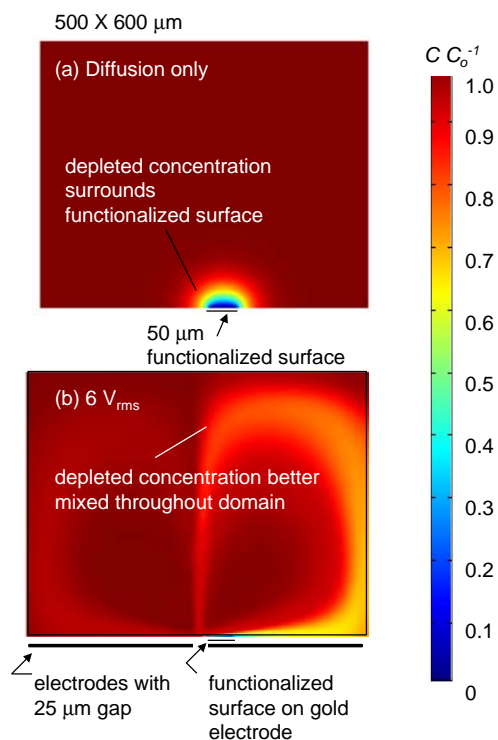

**Fig. 4** Simulation of concentration in a microcavity after 100 sec, subject to binding on the functionalized surface at the wall. In (a) where mass transport is through diffusion only, a depleted cloud (blue) surrounds the functionalized surface, indicating a depletion of analyte. In (b), electrothermal flow circulates the depleted cloud throughout the cavity, and circulates fresh analyte past the functionalized surface, yielding a higher effective concentration close to the functionalized surface.

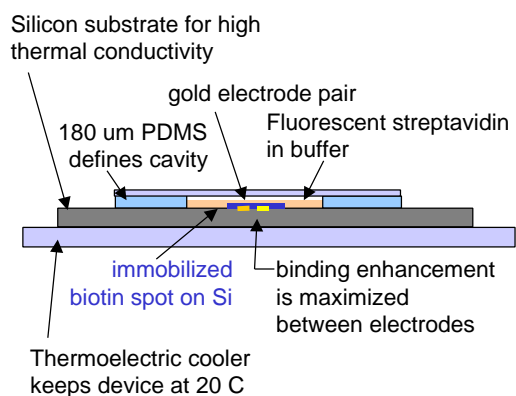

**Fig. 5** Device setup for electrothermally driven microstirring enhancement of streptavidin-biotin heterogeneous assay. A 200 kHz, 10 V<sub>rms</sub> signal is applied across the electrodes to produce electrothermal microstirring. The high thermal conductivity of the silicon base allows a sharp temperature gradient to drive the fluid, while the thermoelectric cooler prevents long-term heating of the system.

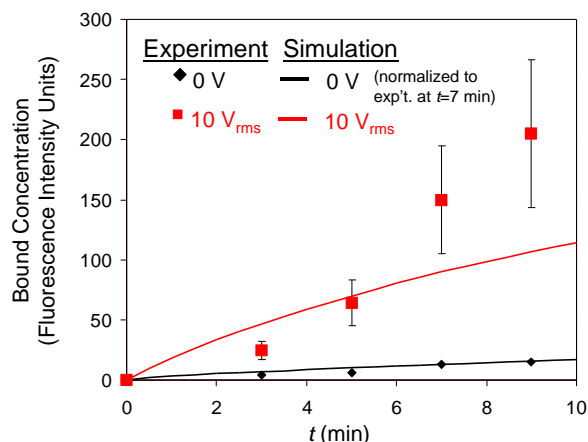

**Fig. 6** Numerical and experimental binding improvement. Bound concentration with time for passive (0V or no stirring; black) and microstirred (10 V<sub>rms</sub>; red). Numerical simulations (solid lines) predict a factor of 7 improvement in the microstirred case, while experimental results show more dramatic improvement factors of 6x (at 6 minutes) to 13x (at 9 minutes). Error bars indicate experimental variation.

the silicon base is maintained at 20 °C with a thermoelectric cooler. This, together with the high thermal conductivity of silicon, provides a “cold” wall, and allows a sufficient temperature gradient to drive the electrothermal flow.

A biotin spot was immobilized on the oxide layer between and around the electrodes. A 150 μm deep microcavity was defined by Parafilm walls and a glass coverslip. The cavity was filled with a fluorescent conjugate of streptavidin and the electrothermal flow was turned on. The reaction was allowed to proceed for 5 minutes, and the slide was rinsed and observed through a fluorescent microscope. The fluorescent intensity is a measure of binding. The flow-enhanced reaction is compared to the zero-flow reaction in Fig. 6, which shows a factor of 6-13 increase in bound concentration with microstirring. Fig. 6 also shows the numerical simulation predicts a factor of 7 increase in bound concentration.

## CONCLUSIONS

Numerical and Experimental results point to a novel method of increasing assay speed by nearly an order of magnitude. This may be particularly useful in portable field-based medical biosensors, or environmental, or food monitoring sensors where quick results are desired. This can be useful for both the warfighter, but also for commercial applications.

## ACKNOWLEDGEMENT

This work is supported by the Institute for Collaborative Biotechnologies through grant DAAD19-03-D-004 from the US Army Research Office.

## REFERENCES

Northrup, M. A., K. F. Jensen, and D. J. Harrison, eds., 2003: Seventh International Conference on Micro Total Analysis Systems. Squaw Valley, CA, Transducers Research Foundation.  
Ramos, A., H. Morgan, N. G. Green, and A. Castellanos, 1998: AC ELECTROKINETICS - A REVIEW OF

FORCES IN MICROELECTRODE STRUCTURES [Review]. *Journal of Physics D-Applied Physics*, **31**, 2338-2353.

Sigurdson, M., D. Z. Wang, and M. C.D., 2005: Electrothermal Stirring for Heterogeneous Immunoassay. *Lab on a Chip*, **5**, 1366-1373.

Stenberg, M. and H. Nygren, 1988: Kinetics of antigen-antibody reactions at solid-liquid interfaces. *Journal of Immunological Methods*, **113**, 3-15.

Wang, D., M. Sigurdson, and C. Meinhart, 2005: Experimental analysis of particle and fluid motion in ac electrokinetics. *Experiments in Fluids*, **38**.